


Positioning Versatility for the Outpatient Setting

The ASC 2000™ Ambulatory
Surgical Care Table



STERIS®

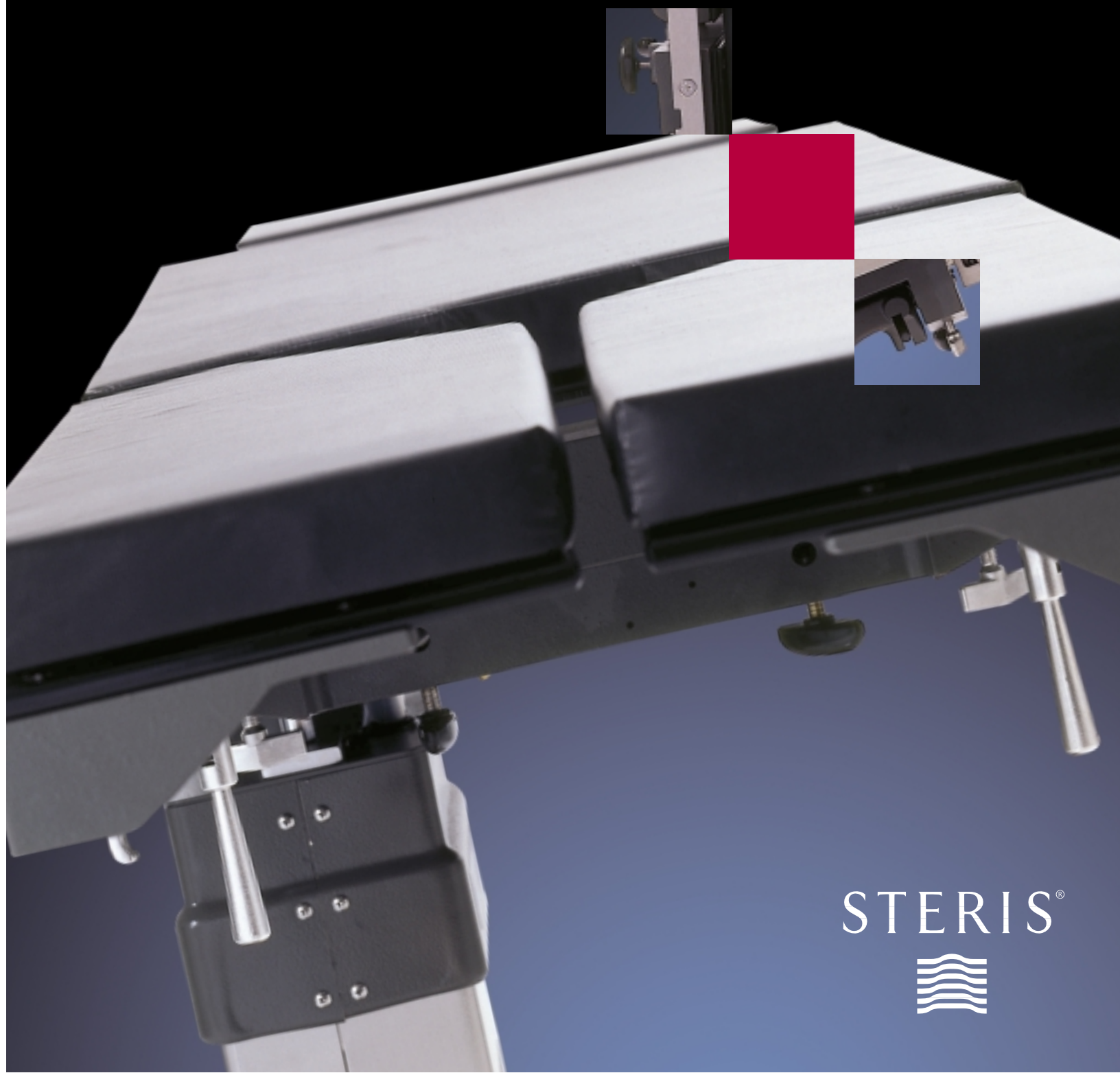

The ASC 2000™ Ambulatory Surgical Care Table

Complete positioning and imaging capability for outpatient surgery


The ASC 2000, a universal, electro-hydraulic operating table, is specially designed to provide complete and flexible positioning for today's outpatient procedures. The unique, offset column design facilitates surgical team access and comfort while it provides superior C-arm access. With its flexibility and ease of use, the ASC 2000 table is designed with the fast paced surgery center in mind.

X-ray Translucence



- The radiolucent tabletop and offset column enable the entire anatomy of patient to be imaged in consistent orientation and without table movement.
- X-ray top is designed so that cassettes can be loaded from the head, foot, or either side for a full range of exposure angles.

Table Movement

- Fully posturing, 6-section tabletop provides flexible positioning for a broad range of procedures.
- Features increased Trendelenburg tilt for greater laparoscopic retraction.
- Battery powered electro-hydraulic actuation with pendant remote hand control for smooth, quiet operation.
- Removable split leg section enables easy leg abduction, mid-line access to abdomen, and perineal approach.

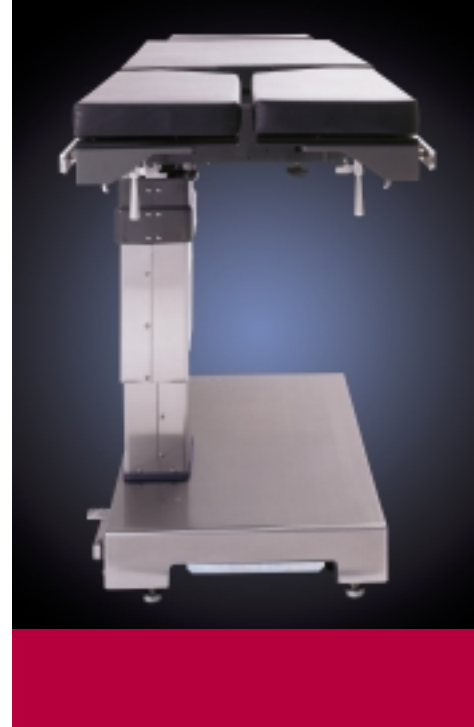


Offset column design provides more than 85% C-arm access, without table movement.

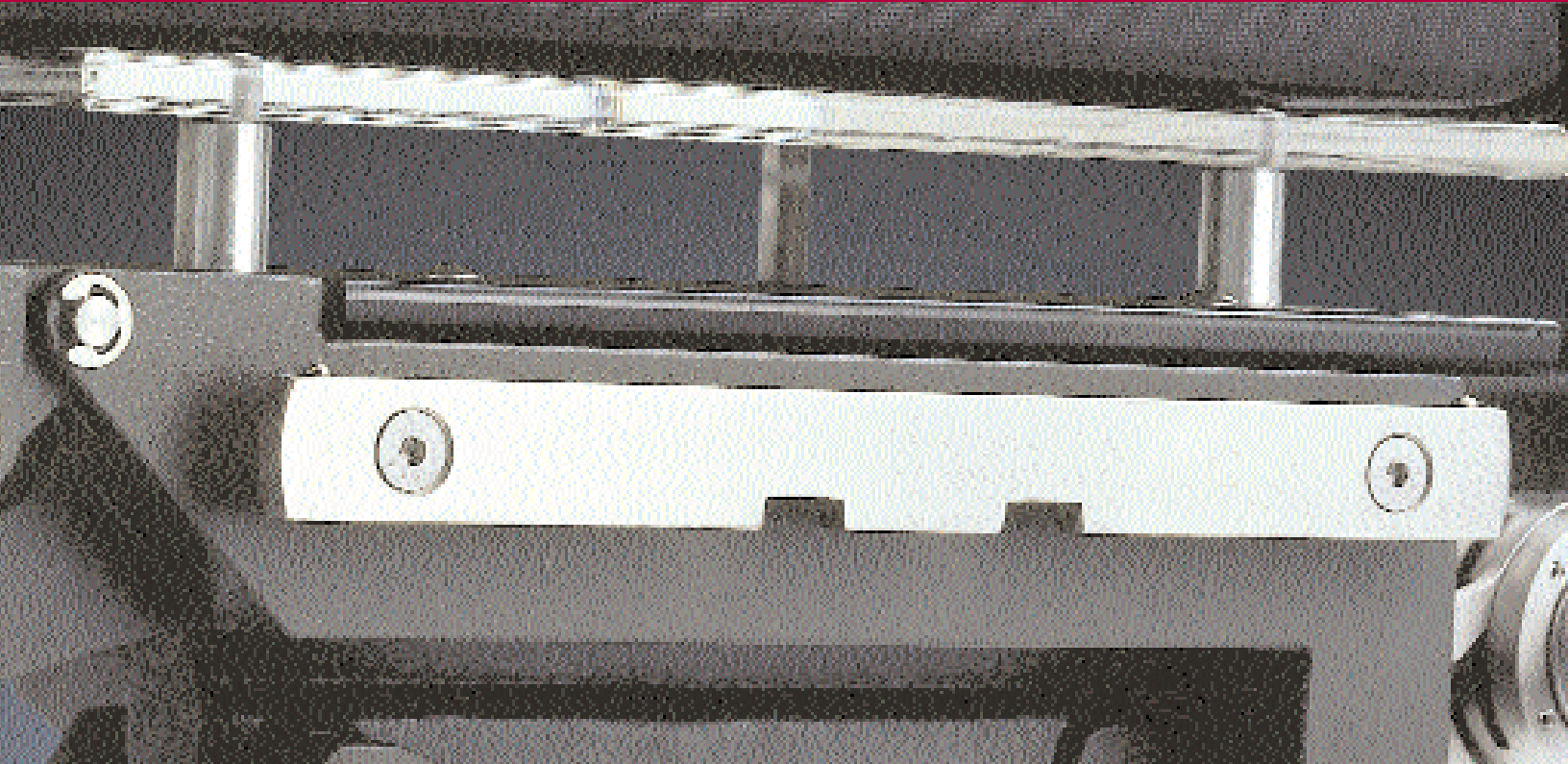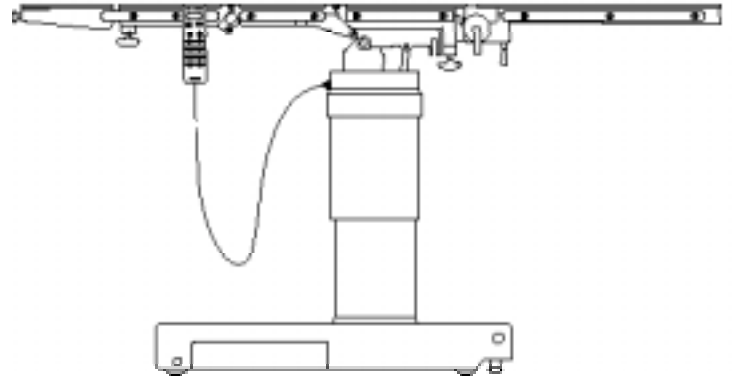


Auxillary control



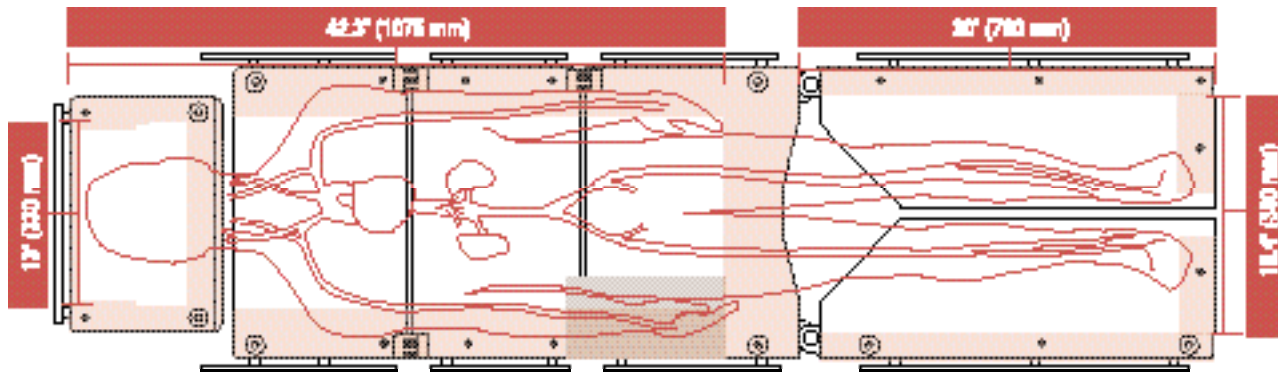

Sized to 6'4" tall patient (97.5th percentile male).

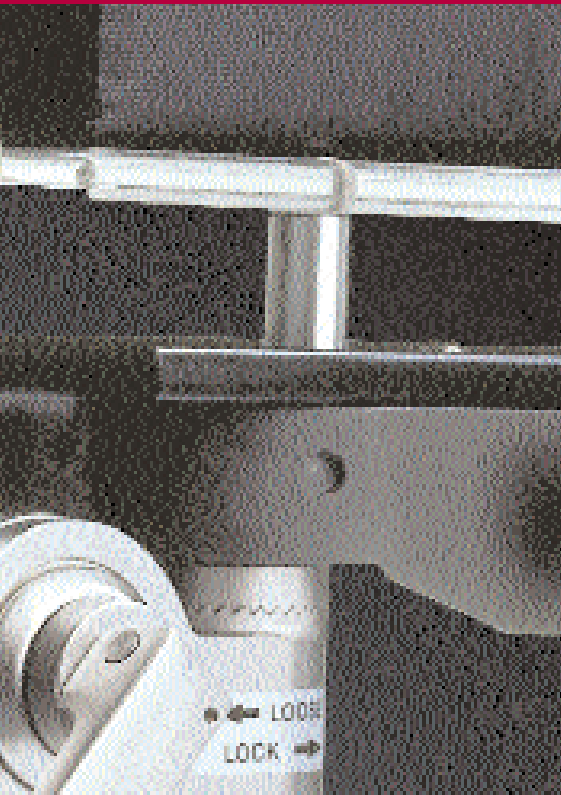

Safety

- The table is durable enough to allow complete articulation of a patient weighing up to 400 lb. (181 kg).
- Self compensating floor locks adjust automatically for uneven surfaces, firmly fixing the table in place.
- Two-inch thick, conductive, pliable, antibiotic table mattress offers patients more safety and comfort, and is easy to maintain.
- Battery power actuation means no cords to clutter the surgical area.

Control

- Convenient remote pendant hand control enables operation from anywhere around the table.
- The 30-second auto standby system helps the surgical team avoid inadvertent operation.
- Auxillary back-up control gives you the ability to articulate the table in the event of hand or foot control failure.
- The optional foot control provides hands-free raising and lowering, Trendelenburg, reverse Trendelenburg, and lateral tilt positioning.

Ease of maintenance

- The beveled edges of the stainless steel cover prevent dirt and fluid from getting into the base.
- The table assembly provides easy service access.
- Available STERIS technical service, including preventative maintenance agreements, helps assure table reliability and long life.

A complete line of accessories supports a broad range of ambulatory procedures



Supine Position

Procedures
Head and Neck
Lymph Node Biopsy
Breast Abnormalities
Hernias
Cholecystectomy
Minilaparotomy
Abdominal

Accessories
X-ray Tops
Armboard (s)
ASC Foot Rests
Conductive Restraint Strap
Ether Screen



Seated Position

Procedures
Breast Abnormalities
Bronchoscopy
Therapeutic Upper Endoscopy
Nissen Fundoplication

Accessories
Multi-Posture Armboard (s)
Clark Sockets or Multi-Accessory
Clamps (Pair)
Conductive Restraint Straps
Hand Control Bracket
Siderail Locks (Pair)
ASC Foot Rests



Lithotomy Position

Procedures
Gynecology
Urology
Endourology

Accessories
Armboard (s)
ASC Drain Pan
Universal Legholders
ASC Pelvic Section Extension
Clark Sockets or Multi-Accessory
Clamps (Pair)
ASC Foot Controller
Ether Screen



Knee/Chest Position

Procedures
Colorectal
Hemorrhoid
Therapeutic Lower Endoscopy

Accessories
Armboard (s)
Morgan Pelvic Rest
Conductive Restraint Strap
Siderail Locks (Pair)
Clark Sockets or Multi-Accessory
Clamps (Pair)
Hand Control Bracket



Modified Fowler/Flex Position

Procedures
Laparoscopic Cholecystectomy

Accessories
Armboard (s)
ASC Foot Rests
Conductive Restraint Strap
Ether Screen



Lateral Flex Position

Procedures
Thoracic

Accessories
Armboard
Multi-Posture Armboard
Lateral Brace Assembly
Conductive Restraint Strap
Siderail Locks (Pair)
ASC Foot Rests



Prone Position

Accessories
Armboard (s)
ASC Foot Rests
Conductive Restraint Strap



Prone Flex Position

Accessories
Armboard (s)
ASC Foot Rests
Conductive Restraint Strap

The ASC 2000™ Ambulatory Surgical Care Table

Technical Data:

Imaging Length of Table	67.5" (1715 mm)
Length:	79.1" (2009 mm)
Width:	19.7" (500 mm)
Width including siderails:	21.5" (547 mm)
Height:	27.6-41.1" (700mm-1045mm)
Trendelenburg/reverse Trendelenburg	45°/20°
Lateral Tilt, L/R	20°/20°
Upper back section Up/Down	40°/55°
Lower back section Up/Down	75°/15°
Legs section Up/Down:	0°/90°
Leg section separate:	0°-90°
Headrest section up/down:	33°/87°
Table Weight:	485 lb (220 kg)

STERIS OFFICES WORLDWIDE

Canada	800 661 3937
France	33 1 55 176000
Germany	49 22 33 69990
Italy	39 02 2130341
Japan	81 78 321 2271
Korea	82 2 517 1517
Latin America	305 442 8202
Singapore	65 68 41 7677
Spain	34 91 658 5920
United Kingdom	44 1256 840400

Publication ID #M2491EN.2004-07, Rev. B
GPSI Printed 07/2004, 5000
©2003-2004 by STERIS Corporation.
All rights reserved. Printed in USA

Technologies to Prevent Infection and Contamination™

STERIS®



STERIS Corporation
5960 Heisley Road
Mentor, OH 44060-1834 ■ USA
440-354-2600 ■ 800-548-4873
www.steris.com